

CT Orbits

Reviewed By: Brian Boldt, DO


Date of Last Revision: August 2026

Date of Last Review: https://docs.google.com/spreadsheets/d/1a1on-GqlksrbCU8zpnnYqisOOebEUGCd_TKeLW5-Clg/edit?gid=0#gid=0

Contact: (866) 761-4200, Option 1

Scan Range: Orbits

Below the orbital ridge through the frontal Sinus

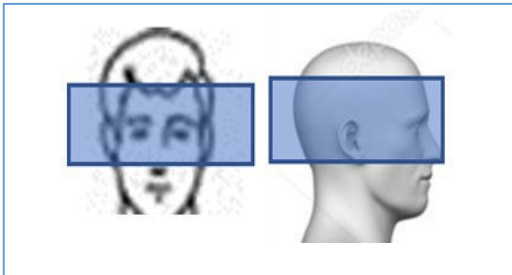

CONTRAST: (Iovue 370) 60 second injection delays.
Injection Volume: 80ml under 250lbs 100ml over 250lbs

Mark sight of symptoms with BB when appropriate.

Reconstructions:

- 1 mm axial soft tissue
- 2 mm coronal soft tissue
- 2 mm sagittal soft tissue
- 1 mm axial bone

*if for trauma add:

- 2 mm coronal bone
- 2 mm sagittal bone

Thresholds:

CTDI:

- 30 mGy 32cm Phantom
- 65 mGy 16cm Phantom